

PRODUCT FACTS

XENON LIGHT SOURCE ILLUMINATORS

- 6000° Kelvin xenon daylight illumination for true tissue rendition
- Low noise fan and highly efficient down-flow cooling system
- 5 year system warranty, 1 year power supply warranty
- Mobile floor stand provides ease of mobility and integral dual hangers secure headlight and fiber optic cable
- Manual variable intensity control allows from 0 – 100% illumination
- Compatible with all IsoLux accessories
- Xenon lamp module provides from 600 to 1,000 hours of lamp life, depending upon usage

300 Watt Xenon

- Compact design 6"H x 11"W x 12"D
- Low-profile, non-conductive, indexing rotary turret designed to accept ACMI/Circon, Olympus, Storz and Wolf type fiber optic cables
- Re-settable digital time meter indicates elapsed number of lamp hours
- UL and CSA approved, CE mark recognition and FDA listed
- Optional storage compartment for loupes, headlights and fiber optic cable

125 Watt Xenon

- Compact design 6"H x 7"W x 11"D
- Outshines any 150 watt halogen or metal halide lamp in the visible spectrum
- Xenon quality at halogen prices
- UL and CSA designed
- Single port design

FIBER OPTIC CABLE

- Fiber optic bifurcated cable, proximal section terminated with Wolf or Universal adapter
- Molded silicone strain-reliever and silicone tubing protects glass fibers
- 2 year fiber optic cable warranty
- Gown clips maintain cable in the desired position

FIBER OPTIC HEADBANDS

Dual Ratchet Headband

- Soft, flexible and ergonomically contoured for maximum comfort with sweat pads of detachable and replaceable plush cushioning foam
- Sculptured ratchet-type knobs provide quick and easy circumferential and crown adjustments prior to or during surgical procedures

Washable Cloth Headband

- Super light weight and softly padded with a secure Velcro closure, ensuring a personalized fit every time
- Hand tool available to remove optics, facilitating the washing of the headband

HEADLIGHT MODULE

- Available with 105° and 90° declination angles to assure the most comfortable surgical stance
- Recessed optical projection lens prevents scratching
- Compact module is designed for use with most loupes
- 10–100 mm variable spot module for precise illumination of the micro surgical field
- 110 mm fixed spot module for large spot illumination
- 5 year headband optics warranty

ORDERING INFORMATION

The management of IsoLux is proud of their products and the greatly improved care surgeons, doctors and dentists around the world are able to provide to their patients because of IsoLux's products. Additionally, management is equally proud of their staff and distributors and the service they are capable of providing.


In the fast paced market place we live in today all too often polite, courteous and efficient service has disappeared. Whether a

customer calls IsoLux directly or a hand-picked distributor that customer will experience service that personalizes their needs and expertly answers their questions.

To place an order for IsoLux equipment, please call 978.774.9136 for your nearest distributor.

IsoLux® llc

Setting the Standard in Medical Fiber Optic Illumination

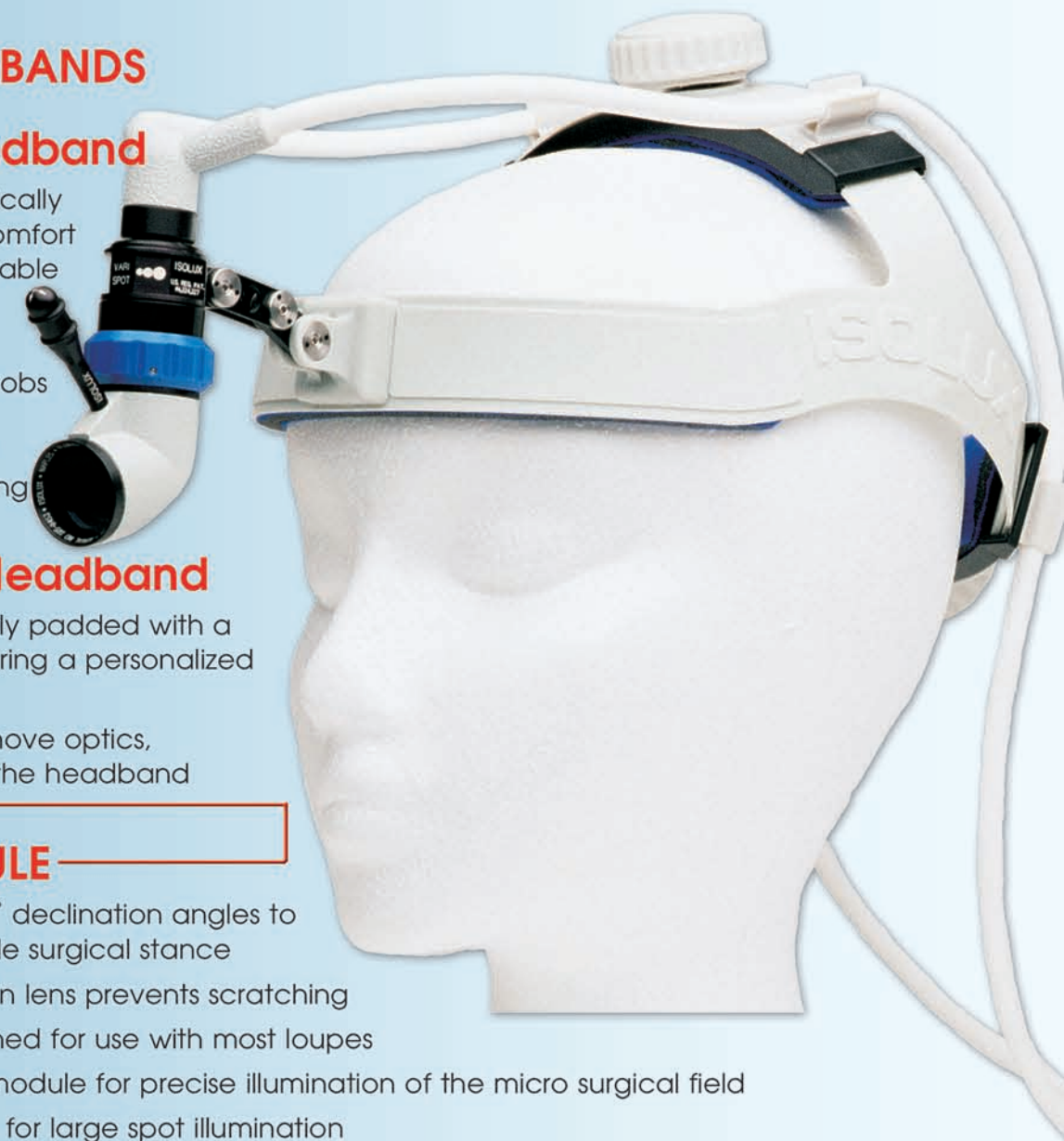

IsoLux® llc